

PAK-6 Antibody

Subcategory: Rabbit Polyclonal Antibody

Cat. No.: 251927

Unit: 0.1 mg

Description:

Serine/threonine-protein kinase PAK 6 (PAK-6) is a protein kinase that acts on a variety of targets. PAK-6 interacts tightly with GTP-bound but not GDP-bound CDC42/p21 and RAC1, and also interacts with the androgen receptor. PAK-6 is autophosphorylated when activated by CDC42/p21.

Isotype: Rabbit Ig

Applications: E, WB

Species Reactivity: H


Format: Each vial contains 0.1 mg IgG in 0.1 ml (1 mg/ml) of PBS pH7.4 with 0.09% sodium azide. Antibody was purified by Protein-G affinity chromatography.

Alternate Names: Serine/threonine-protein kinase PAK 6; p21-activated kinase 6; PAK-6; PAK-5; PAK6; PAK5

Accession No.: Q9NQU5

Antigen: KLH-conjugated synthetic peptide encompassing a sequence within the N-term region of human PAK-6.

Application Notes: E: 1:1,000-1:20,000; WB: 1:100-1:500

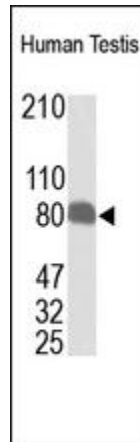

The PAK-6 Antibody is used in Western blot to detect PAK-6 in human fetal testis lysates.

Storage: Store at -20°C. Minimize freeze-thaw cycles. Product is guaranteed one year from the date of shipment.

For research use only, not for diagnostic or therapeutic procedures.