

TNFRSF18 Antibody

Subcategory: Rabbit Polyclonal Antibody

Cat. No.: 251542

Unit: 0.1 mg

Description:

Tumor necrosis factor receptor superfamily member 18 (TNFRSF18) is a receptor for TNFRSF18. TNFRSF18 is involved in interactions between activated T-lymphocytes and endothelial cells and in the regulation of T-cell receptor-mediated cell death. TNFRSF18 mediates NF-kappa-B activation via the TRAF2/NIK pathway. TNFRSF18 is expressed in lymph node, peripheral blood leukocytes and weakly in spleen.

Isotype: Rabbit Ig

Applications: E, IHC, WB

Species Reactivity: H, M, R


Format: Each vial contains 0.1 mg IgG in 0.1 ml (1 mg/ml) of PBS pH7.4, 25% glycerol, 0.1 mg/ml BSA with 0.09% sodium azide. Antibody was purified by Protein-A affinity chromatography.

Alternate Names: Tumor necrosis factor receptor superfamily member 18; Glucocorticoid-induced TNFR-related protein; Activation-inducibleTNFR family receptor; TNFRSF18; AITR; GITR

Accession No.: Q9Y5U5

Antigen: KLH-conjugated synthetic peptide encompassing a sequence within the C-term region of human TNFRSF18 .

Application Notes: E: 1:500-1:1,000; WB: 1:100-1:500; IHC: 1:100-1:500

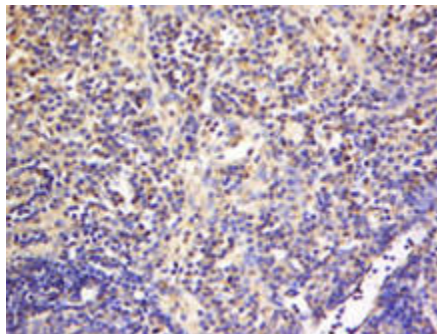

TNFRSF18 staining in mouse lymph node. Formalin-fixed paraffin-embedded mouse lymph node tissue is stained with TNFRSF18 Antibody (Cat. No. 251542) used at 1:200 dilution.

Storage: Store at -20°C. Minimize freeze-thaw cycles. Product is guaranteed one year from the date of shipment.

For research use only, not for diagnostic or therapeutic procedures.