

SMPDL3B Antibody

Subcategory: Rabbit Polyclonal Antibody

Cat. No.: 251962

Unit: 0.1 mg

Description:

Acid sphingomyelinase-like phosphodiesterase 3b (SMPDL3B) is a member of a multi-gene family and shares motifs with a large family of metallo-phosphoesterases.

Isotype: Rabbit Ig

Applications: E, IHC, WB

Species Reactivity: H


Format: Each vial contains 0.1 mg IgG in 0.1 ml (1 mg/ml) of PBS pH7.4, 1% BSA with 0.09% sodium azide. Antibody was purified by Protein-G affinity chromatography.

Alternate Names: Acid sphingomyelinase-like phosphodiesterase 3b; ASM-like phosphodiesterase 3b; SMPDL3B; ASML3B

Accession No.: Q92485

Antigen: KLH-conjugated synthetic peptide encompassing a sequence within the center region of human SMPDL3B.

Application Notes: E: 1:10,000-1:40,000; WB: 1:100-1:500

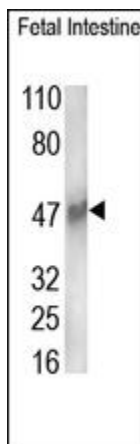

The SMPDL3B Antibody (Cat. No. 251962) is used in Western blot to detect SMPDL3B in human fetal intestine lysates.

Storage: Store at -20°C. Minimize freeze-thaw cycles. Product is guaranteed one year from the date of shipment.

For research use only, not for diagnostic or therapeutic procedures.